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INSTRUCTION MANUAL

Read this manual carefully before using the device and save it for future reference. To reduce possible personal injury, keep in mind the following instructions:

CAUTION

- ▶ Do not charging it in bath room.
- ▶ Do not break off products or re-change battery, if need change battery,ask for help from stuffs.
- ▶ Do not change USB cable line,send it to after-service department if need.
- ▶ Checking if the voltage of product adapt to your local voltage before recharging it.
- ▶ One set includes handle, cable type C and toothbrush heads. Any part breaks down,replace it to avoid dangers. (only replace it by assorted replacements from us)Giving it away when it can not replace.
- ▶ Put product out of children&infant' reach.Do not play it as a play toy.
- ▶ Please ask your doctor before using it if you have a operation in 2 months.
- ▶ Do not using it for other purpose except for oral health care
- ▶ Clearing up brush heads after using when your toothpaste with hydrogen peroxide and dicarbonate.
- ▶ Do not using it when taking shower.

Turning it off when not using it

CHARGING BATTERY

- ▶ Running out all the energy for the first time to use. Recharging it for next time. Running out all energy for every month to protect the battery life.
- ▶ Recharging for 20-24h after 6-months no use.
- ▶ Long time no recharging, the indicator light will on after about 6-10min.

LED	COLOR
Working	Display in white light
Low power protection	Red light flash
Low power	Red light twinkling
Charging	White light flash
Full battery	All led light in white

FUNCTION DESCRIPTION

- ▶ Modes: press the switch for 8s to change into 3 different modes; *sensible, blanquear, limpiar, apagado.*
- ▶ OFF/ON:Press the switch for twice starting toothbrush.Press switch for 8s the toothbrush turn off.
- ▶ Pause when start:after start toothbrush 3s,the frequency will come to normal.
- ▶ Power check:press switch shortly to check power.

MODE

Sensitive-white-clean, In any modes, it will pause for 0.3s to remind you to brush next area. The total brushing time will be 2 mir.

COMPONENTS

1. Brush head
2. Switch
3. Function *Sensitive*
4. Function *White*
5. Function *Clean*
6. USB input (for cable type C)


BODY	
Voltage and power	DC 3,7V 1W
Working time	45 days approx.
Frequency	34200 - 40200RPM/min
Battery Type	Lithium battery
Charging time	3h

PACKING LIST

- Body: x1
- Cable Type-Cline: x1
- Brush head: x2
- User manual: x1

ERROR DIAGNOSIS

ERROR	REASONS	ACTION
Frequency Weaken	• Too hard to press toothbrush on teeth when using • Toothbrush heads fix too firmly on toothbrush body.	• 1 mm gap between toothbrush heads and toothbrush body is necessary. • Put toothbrush on teeth lightly.
Full charge but short working time	• Low power • Not connect correct with cable line.	Connect cable Type C correct with bottom charge jack, and check
Toothbrush not workable	• Charging time is not enough. • Not using it for over 6 months.	• Re-charge toothbrush • Buy new battery from after-service department and replace it.

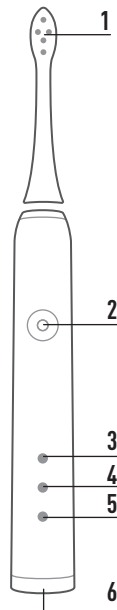MAINTENANCE

Clearing brush head after using and drying it by air.

Clear body by soft cloth
Do not put it in the water
Do not clear it by 40 C above water or acid solvent

Replace brush head
Replace brush head for 3 month as dentist suggests

TIMING

- Brushing tooth for 2min as dentist suggest.
- The device will remind user to change brushing area in every 30 seconds.
- The device will record the mode after every single use.

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ENGLISH

WASTE

This symbol on the product or in the instructions means that your electrical and electronic equipment must be disposed at the end of its useful life in a specialized container; Do not dispose it in the usual household waste container. In the EU there is special waste collection systems for subsequent recycling. For more information, contact with the local authority or the retailer from whom you purchased the product.

WARRANTY

ELEKTRO-3, S.C.C.L. certifies that its products comply with the quality requirements and HAVE THE LEGAL GUARANTEE according to the regulations in force in each country. In case of observing any defect or problem in the product, please contact our Technical Assistance Service (SAT) by calling +34 610 699 257.

THE INVOICE IS THE PRODUCT WARRANTY

THE MANUFACTURER DECLINES ALL RESPONSIBILITY DUE TO AN INCORRECT INSTALLATION OF THE APPLIANCE.

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